


Polypropylene Ductless Fume Hood

Topair's Polypropylene Ductless Fume Hoods provide an ultimate protective work environment for lab staff working with acids and harsh chemicals.

The electrical and mechanical components are manufactured by leading global companies, such as EBM Germany and AAF USA. The products comply with relevant international standards.

TopAir Polypropylene Ductless Fume Hoods are customized to the requirements of each client.

- Welded white polypropylene structure
- Built-in sealed polypropylene worktop
- Easily dissembled back wall
- Tempered glass sliding front window
- Monitor indicating fan's total working time, for tracking and filter replacement purposes
- Top quality, quiet fan manufactured by EBM Germany
- Eco-friendly, cost-effective 600-800 LUX LED lighting separate from the fuming hall.
- Wind speed at $0.6+0.1$ m/s, $120\pm$ FPM
- Easy filter replacement
- Top filtration unit including carbon filter
- User-friendly digital control system including fan speed control, and displays for air pressure, signal light, operation time and alarms.

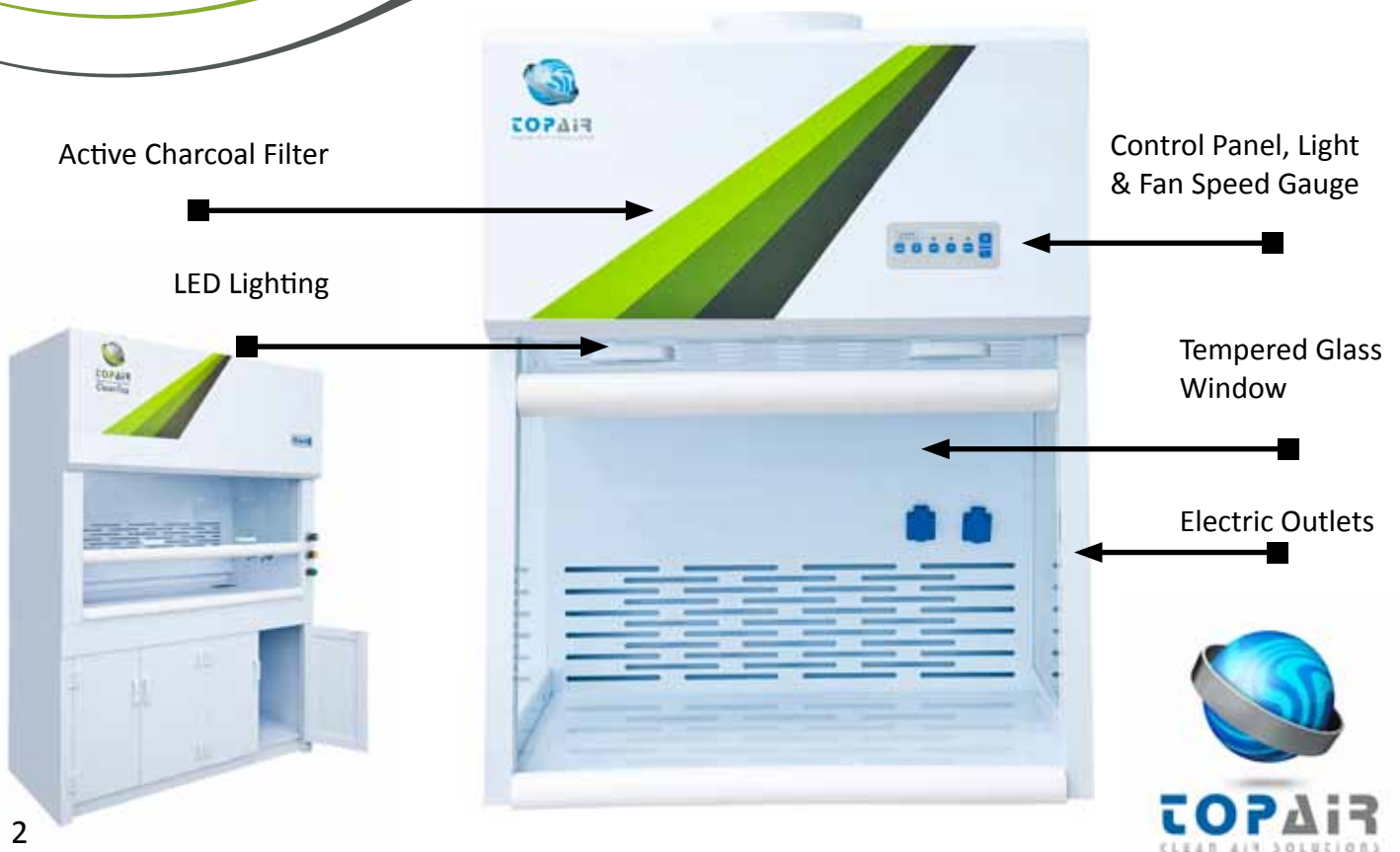

Models

Spec/ Model	CF-060-A-PP	CF-090-A-PP	CF-120-A-PP	CF-160-A-PP	CF-180-A-PP
Outer Dimensions	600 x 700 x 1050 mm 23.62 x 27.5 x 41.3"	900 x 700 x 1050 mm 35.4 x 27.5 x 41.3"	1200 x 700 x 1050 mm 47.24 x 27.5 x 41.3"	1600 x 700 x 1050 mm 62.99 x 27.5 x 41.3"	1800 x 700 x 1050 mm 70.67 x 27.5 x 41.3"
Workspace (W x D x H)	540 x 640 x 660 mm 21.26 x 25.2 x 26"	840 x 640 x 660 mm 33.07 x 25.2 x 26"	1140 x 640 x 660 mm 44.88 x 25.2 x 26"	1540 x 640 x 660 mm 60.63 x 25.2 x 26"	1740 x 640 x 660 mm 66.50 x 25.2 x 26"
Production / Test Standard	EN-14175 / CE / ASHRAE 110-1995 / ISO-9001:2008				
Air Velocity	0.6±0.1m/s, 120±20 FP				
Cabinet Material	Welded white polypropylene structure with built-in sealed polypropylene worktop				
Noise Level	<52dB	<52dB	<54dB	<60dB	<62dB
	(Tested 20 cm from the work table, 1.2m above ground)				
Power Supply Options	110 / 220V 50/60 Hz, Single phase				
Illumination	>600-800 LUX, Eco-friendly LED lighting				
Filter	Charcoal Filter/ multi-gas filter/HEPA				

Accessories

Spec/ Model	CF-060-A-PP	CF-090-A-PP	CF-120-A-PP	CF-160A-PP	CF-180-A-PP
Stand WxDxH	CF-060-ST 600 x 700 x 850 mm 23.6 x 27.56 x 33.46"	CF-090-ST 900 x 700 x 850 mm 35.43 x 27.56 x 33.46"	CF-120-ST 1200 x 700 x 850 mm 47.24 x 27.56 x 33.46"	CF-160-ST 1600 x 700 x 850 mm 63 x 27.56 x 33.46"	CF-180-ST 1800 x 700 x 850 mm 70.86 x 27.56 x 33.46 "
UV Light	CF-060-UV	CF-090-UV	CF-120-UV	CF-160-UV	CF-180-UV
Separate Table	CF-060-UB	CF-090-UB	CF-120-UB	CF-160-UB	CF-180-UB

Web: www.topairsystems.com Email: sales@topairsystems.com

Office - USA: 3182 Monterey Dr., Merrick, NY 11566, USA

Tel : +1-718-841-8886 Email: sales@topairsystems.com

**Office - UK: Evolution Testing & Analytical services (UK) Ltd,
Elstree House, Elstree Way, Borehamwood, Herts WD6 1SD, UK**

Tel:+44-203-1374012 Email: sales@topairsystems.com

Copyright © 2014 TopAir Systems. All Rights Reserved.