

# Educational Ductless Fume Hood with All-Round Clear Glass


Topair's Polypropylene Ductless Fume Hoods provide a safe work environment for lab staff working with acids and harsh chemicals.

The electrical and mechanical components are manufactured by leading global companies, such as AAF USA. The products are **EN-14175 / CE / ASHRAE 110-1995 certified**.


TopAir's Polypropylene Ductless Fume Hoods are customized to the requirements of each client.

**Optional Features**

PPS V2      ESD Protected



- 360° transparency for exceptional visibility
- Polypropylene structure, high chemical resistance
- Built-in sealed polypropylene worktop
- Optional stand
- Easily dissembled back wall
- Tempered glass sliding front window
- Monitor displays fan's total operation time, for tracking and filter replacement purposes
- Top quality quiet fan
- Eco-friendly, cost-effective 800 LUX LED lighting separated from the work area
- Wind speed at 0.5±0.1 m/s, 100±20 FPM
- Convenient front access for filter replacement
- Top filtration unit - variety of HEPA & carbon filters
- User-friendly digital control system
- EN-14175 / CE / ASHRAE 110-1995 certified

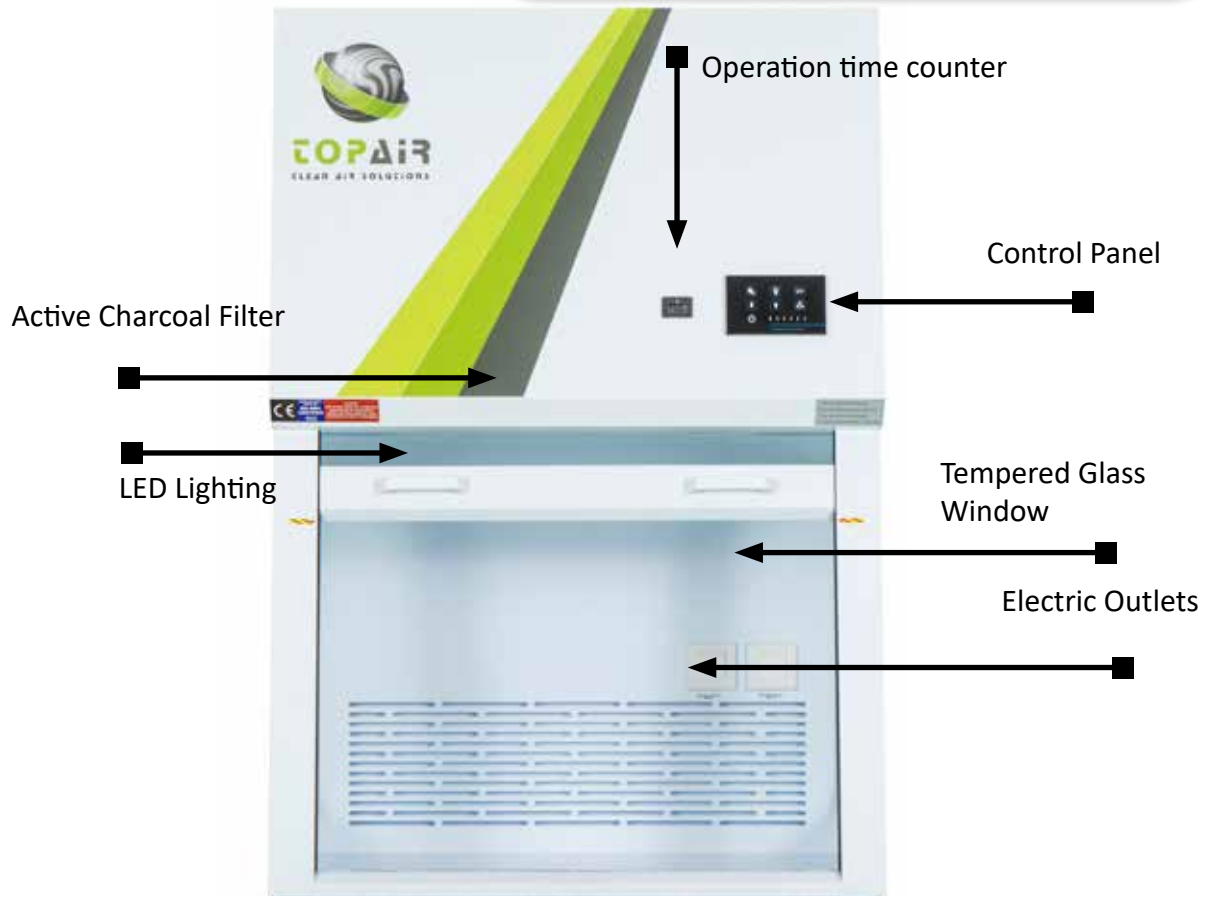

# Models

**CUSTOM SIZES AVAILABLE!**

Spec/ Model	CF-060-CB	CF-090-CB	CF-120-CB	CF-150-CB	CF-180-CB
Outer Dimensions	600 x 750 x 1223 mm	900 x 750 x 1223 mm	1200 x 750 x 1223 mm	1500 x 750 x 1223 mm	1800 x 750 x 1223 mm
W x D x H	23.62 x 29.5 x 48"	35.4 x 29.5 x 48"	47.24 x 29.5 x 48"	59 x 29.5 x 48"	70.8 x 29.5 x 48"
Workspace (W x D x H)	585 x 610 x 695 mm	885 x 610 x 695 mm	1185 x 610 x 695 mm	1485 x 610 x 695 mm	1785 x 610 x 695 mm
	23 x 24 x 27.3"	34.8 x 24 x 27.3"	46.6 x 24 x 27.3"	58.4 x 24 x 27.3"	70.2 x 24 x 27.3"
Front Sash Max. Opening	570 mm / 22.4"				
Production / Test Standard	EN-14175 / CE / ASHRAE 110-1995				
Air Velocity	0.5±0.1 m/s, 100±20 FPM				
Hood Material	Welded white polypropylene structure with built-in sealed polypropylene worktop				
Noise Level	<52dB	<52dB	<54dB	<60dB	<62dB
	(Tested 20 cm from the work table, 1.2m above ground)				
Power Supply	110 / 220V, 50/60 Hz, Single phase				
Illumination	800 LUX, Eco-friendly LED lighting				
Filter	Charcoal Filter/ multi-gas filter/HEPA				

## Accessories

Description	Model
Metal stand	CF-size-ST
Polypropylene base cabinet	CF-size-BS
UV light	CF-UV
Gas tap	CF-GTAP
Water tap	CF-WTAP
Polypropylene cup sink	CF-PP-SINK
Polypropylene sink 30 x 40	CF-SINK-3040
Power outlet installed	CF-SOCKET

Tel : +1-855-6-TOPAIR (+1-855-686-7247) Email: [sales@topairsystems.com](mailto:sales@topairsystems.com)

Web: [www.topairsystems.com](http://www.topairsystems.com) | Headquarters USA: 8912 68th Avenue Forest Hills New York 11375, USA