

Polypropylene Biosafety Cabinet

TopAir's Class II Biological Safety Cabinets protect staff, the environment, and sensitive work processes in which biological agents are applied, typically in the biological and microbiological industries.


The product offers a high level of contamination protection, based on two advanced HEPA filters operating at a typical efficiency of 99.9995% @ 0.3 um.

The polypropylene structure offers an optimal solution as a solid, easily-cleaned, high resistance material.

The cabinet is equipped with a smart, safe and elegant touch-screen control system that protects the worker and provides alerts on periodic maintenance actions and device replacement.

All components are low energy consumption, with LED lighting and EC fan motor. The system also has programmable "green" night mode, that shuts off all unnecessary electricity consumers and sets the vital components at the required safety level.

- Polypropylene structure
- 6 mm tempered glass side walls
- 304 stainless steel work surface and spill tray
- Two H14 HEPA filters
- Advanced EC fan with 304 stainless steel housing
- 55 dba noise level
- ISO 5/ CLASS 100 cleanliness level
- Smart, programmable advanced touch screen control system
- Technician calibration screen
- Maintenance alarms
- Technical faults alarms
- Timers and counters management screen
- Germicidal water proof UV light system and safety interlock mechanism
- 6 mm double layer safety front glass with electrical motion system
- Programmable economical night mode
- Economical LED light

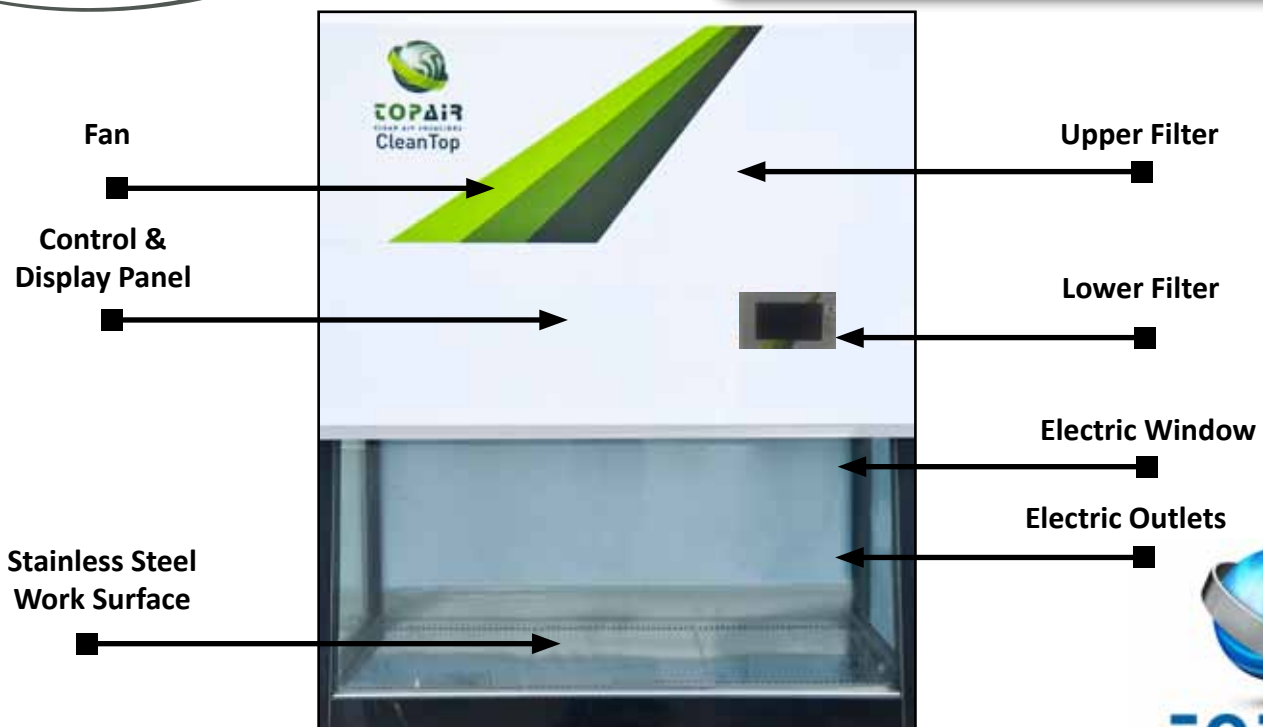

Models

Spec/ Model	BO-2A-090PP	BO-2A-120PP	BO-2A-150PP	BO-2A-180PP
Outer Dimensions W x D x H	900 x 800 x 1300 mm 35.4x 31.5 x 51.2"	1200 x 800 x 1300 mm 47.2 x 31.5 x 51.2"	1500 x 800 x 1300 mm 59 x 31.5 x 51.2"	1800 x 600 x 1300 mm 70.9 x 31.5 x 51.2"
Workspace (W x D x H)	850 x 650 x 550 mm 33.46 x 25.6 x 21.65"	1150 x 650 x 550 mm 45.27 x 25.6 x 21.65"	1450 x 650 x 550 mm 57 x 25.6 x 21.65"	1750 x 650 x 550 mm 68.9 x 25.6 x 21.65"
Production/ Test Standard	CE / ISO-9001:2008/EN12469			
Downflow Velocity	0.45 m/s, 90 Fpm			
Inflow velocity	0.5 m/s, 100 fpm			
Cleanliness level	Class 100/ISO 5			
Cabinet Material	Welded white polypropylene structure with stainless still 304 worktop			
Noise Level	<52dB	<52dB	<54dB	<60dB
	(Tested 20 cm from worktable, 1.2m above ground)			
	115 / 230V 50/60 Hz, Single phase			
Illumination	>600-800 LUX, Eco-friendly LED lighting			
Filters	HEPA/ULPA			

Web: www.topairsystems.com **Email:** sales@topairsystems.com

Office - USA: 3182 Monterey Dr., Merrick, NY 11566, USA

Tel : +1-718-841-8886 **Email:** sales@topairsystems.com

Office - UK: Evolution Testing & Analytical services (UK) Ltd,

Elstree House, Elstree Way, Borehamwood, Herts WD6 1SD, UK

Tel:+44-203-1374012 **Email:** sales@topairsystems.com

Copyright © 2014 TopAir Systems. All Rights Reserved.