Ductless Fume Hood-PRO



NEW

Topair's Polypropylene Ductless Fume Hood - PRO provides a safe work environment for lab staff working with acids and harsh chemicals.

The hood includes an advanced VAV (Variable Air Volume) system. The VAV system measures the product's air velocity using a high quality sensor, and adjusts the air velocity speed to the relevant standard.


Maximal energy savings are enabled by flexibly adjusting the fan speed (high/low) as needed. The system features a high safety level, displaying real time air velocity and providing alarms for low velocity levels, thus reducing the fume hood's noise level.

Advanced Operation System

- Color 9" touch screen
- Hour counter for filter
- Filter replacement alarm for HEPA/Carbon filters
- Sensor for temperature and humidity

The products are EN-14175 / CE / ASHRAE 110-1995 certified and comply with AFNOR NFX 15-211 (Class 1 & 2)

- Automatic airflow control (VAV) with 9" color touch screen, various operation speeds, configurable alarms, visual & audio alarms.
- Chemical sensor alarm
- Polypropylene structure, high chemical resistance
- Built-in sealed polypropylene worktop or choice of epoxy, stainless steel, ceramic, Trespa
- Optional stand,
- High efficiency EC fan
- Tempered glass sliding front window
- 800 LUX LED lighting
- Air velocity 0.5 m/s, 100 FPM
- Convenient front access for filter replacement
- Variety of HEPA & carbon filters
- Exhaust point for connecting a 2nd filter or an external exhaust.
- XL model with 2 internal filters is also available
- Complies with AFNOR NFX 15-211 (Class 1 & 2)
- EN-14175 / CE / ASHRAE 110-1995 certified

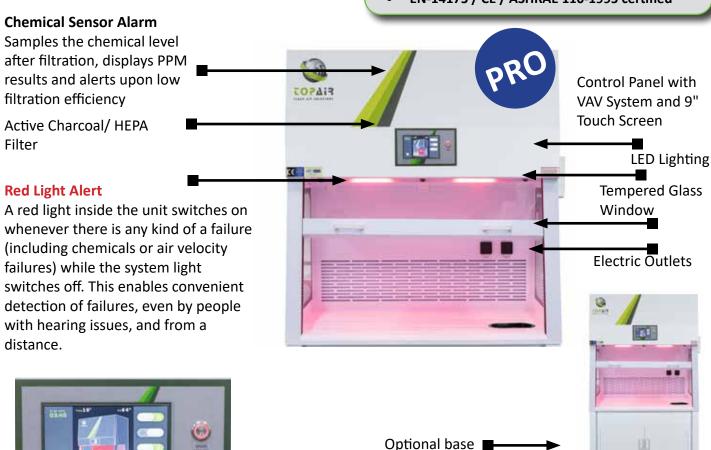

Models

Spec/ Model	CF-060-PRO	CF-090-PRO	CF-120-PRO	CF-150-PRO	CF-180-PRO
Outer	600 x 750 x 1223	900 x 750 x 1223	1200 x 750 x 1223	1500 x 750 x 1223	1800 x 750 x 1223
Dimensions	mm	mm	mm	mm	mm
WxDxH	23.62 x 29.5 x 48"	35.4 x 29.5 x 48"	47.24 x 29.5 x 48"	59 x 29.5 x 48"	70.8 x 29.5 x 48"
Dimensions for XL model with 2 internal filters W x D x H	600 x 750 x 1500 mm	900 x 750 x 1500 mm	1200 x 750 x 1500 mm	1500 x 750 x 1500 mm	1800 x 750 x 1500 mm
For XL models add –XL to P/N	23.62 x 29.5 x 59"	35.4 x 29.5 x 59"	47.24 x 29.5 x 59"	59 x 29.5 x 59"	70.8 x 29.5 x 59"
WE CAN CUSTOMIZE TO ANY SIZE - EVEN A SINGLE UNIT! CONTACT US FOR DETAILS					
Workspace	585 x 590 x 695 mm	885 x 590 x 695 mm	1185 x 590 x 695 mm	1485 x 590 x 695 mm	1785 x 590 x 695 mm
(W x D x H)	23 x 23.2 x 27.3"	34.8 x 23.2 x 27.3"	46.6 x 23.2 x 27.3"	58.4 x 23.2 x 27.3"	70.2 x 23.2 x 27.3"
Front Sash Max. Opening	495 mm / 19.5"				
Test Standard	EN-14175 / CE / ASHRAE 110-1995				
Air Velocity	0.5±0.1 m/s, 100±20 FPM				
Hood Material	Welded white polypropylene structure with built-in sealed polypropylene worktop				
Noise	<52dB	<52dB	<54dB	<60dB	<62dB
(Tested 20 cm from the work table, 1.2m above ground)					
Power Supply	110 / 220V, 50/60 Hz, Single phase				
Illumination	800 LUX, Eco-friendly LED lighting				
Filter	Carbon filters for solvents, acids, bases, formaldehyde, multi gas (blend of materials) /HEPA/ULPA				

Accessories

Description	Model		
Metal stand	CF-size-ST		
Polypropylene base cabinet	CF-size-BS		
UV light	CF-size-UV		
Gas tap	CF-GTAP		
Water tap	CF-WTAP		
Polypropylene cup sink	CF-PP-SINK		
Polypropylene sink 30 x 40	CF-SINK-3040		
Power outlet installed	CF-SOCKET		

Tel: 1-855-6-TOPAIR International: +1-855-686-7247 Email: sales@topairsystems.com www.topairsystems.com | Headquarters: 300 First Avenue, Suite 102, Needham, MA 02494 USA